
ACTUATION SYSTEMS for PATIENT LIFTING APPLICATIONS

TiMOTION Provides Flexibility, Customized Products, and Tailor-Made Services



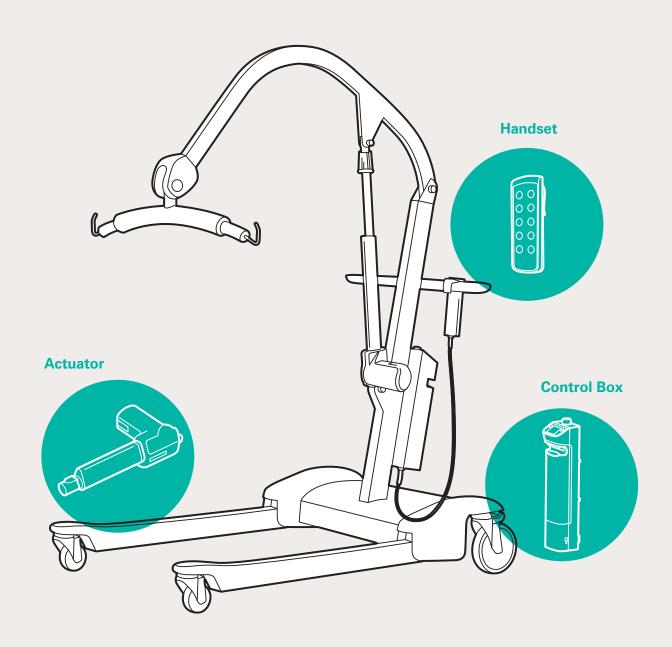

AN INTELLIGENT ACTUATION SYSTEM FOR PATIENT LIFTS

TiMOTION created an optimum solution to meet specific patient needs and to secure a safe and comfortable transfer.

The TiMOTION actuation system also provides a safe and efficient working environment for the caregivers.







AN INTELLIGENT CONTROL SYSTEM TC12 SERIES

- · IP protection: IPX4
- · EN 10535: 2006 and IEC60601-1:2005 compliant
- · Multiple options of display: without / black-white / color
- · Includes battery box
- · With wall mount battery charger or battery charger dock
- · Regular maintenance / overload alarm notification
- · Battery alarm when battery with lower capacity remaining
- · Easy installation

Key Features of the UI on TC12*


Note: here's an example of the TC12 with color-display


Drawing

Standard Dimension (mm)

Battery type




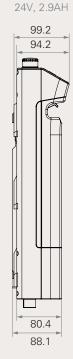


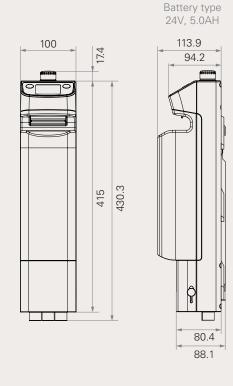


Operation Status

Control open / close, and up / down movements

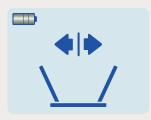

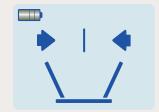



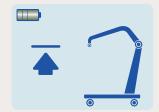

Notification

Alarm for low battery capacity, overloading, and safety switch to prevent over extending









Battery Status










Note: for more detailed user guide information, please contact TiMOTION sales team.



Solution Overview

Control box: TC12 Handset: TH10

Actuators for lifting: TA1 or TA12

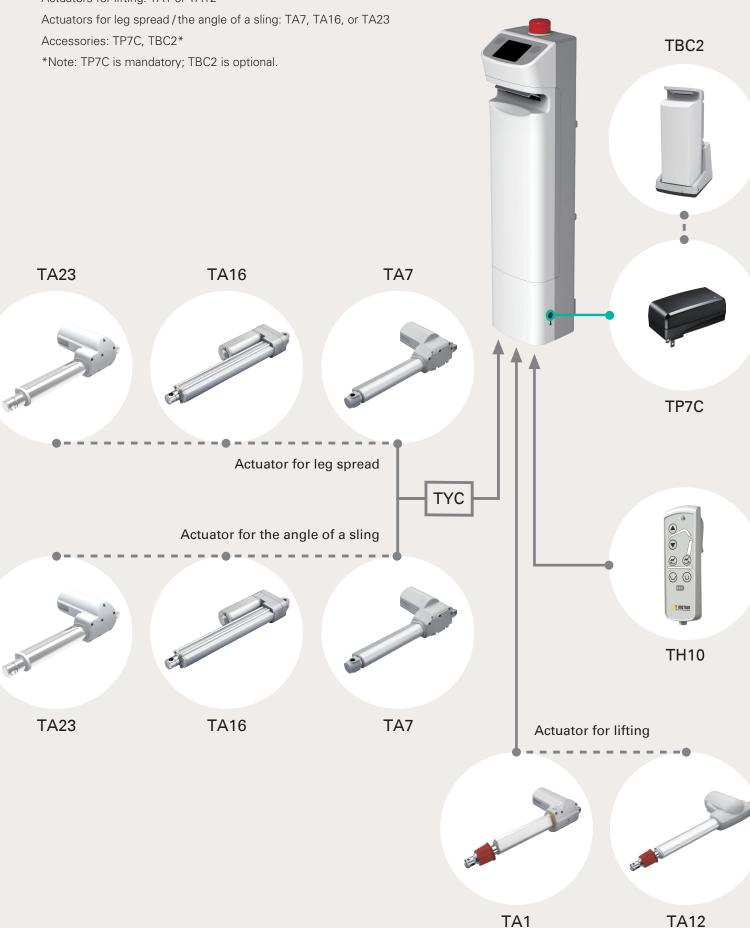



ACTUATOR for LIFTING TA1 SERIES

· Load: 6,000 / 8,000N in push

· Protection class: IPX4

· Certificate: IEC 60601-1 compliant




BATTERY CHARGER for TC12 TP7C SERIES

· The TP7C must always be ordered separately when ordering TC12

HANDSET

TH10 SERIES

- · Indicator lights show battery status
- · Controls lifting and lowering
- · Controls lower leg movement
- · Alarm notifications for low battery capacity and overloading

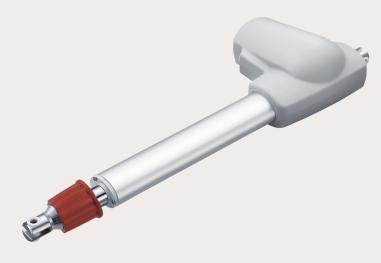

ACTUATOR for LIFTING TA12 SERIES

· Load: 10,000 / 12,000N in push

· Protection class: IPX4

· Certificate: IEC 60601-1 compliant



CHARGING CRADLE for TC12 TBC2 SERIES

· TBC2 servers as a battery charger (optional)



· For single motor




· For two motors



· For three motors

Installation Guide


1. Install the bracket on the hoist frame

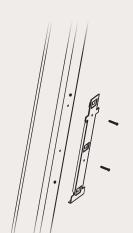


2. Install the TC12 on the bracket




3. Make sure all fasteners are tightened properly

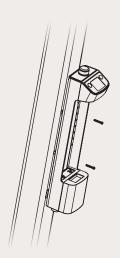

4. Install the battery

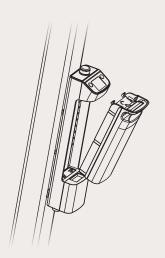

5. Installation is complete

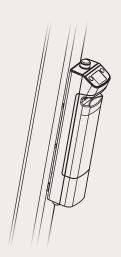

Headquarter

TiMOTION Technology Co., Ltd.
Tel 886 2 2219 6633
Email sales@timotion.tw

TiMOTION Europe

Tel 33 (0)1 74 82 50 51 Email sales@timotion.fr

TiMOTION Kunshan

Tel 86 512 5526 0735 Email sales@timotion.tw

TiMOTION Japan

Tel 81(0)6 6763 1110 Email sales@timotion.jp

TiMOTION Korea

Tel 82 31 745 1060 Email sales@timotion.co.kr

TIMOTION USA

Tel (704) 708 6924 Email sales@timotionusa.com

